

Purification of His-tag proteins

User manual

Protino[®] Ni-NTA Agarose Protino[®] Ni-NTA Columns 1 mL Protino[®] Ni-NTA Columns 5 mL

October 2013/Rev. 04



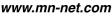


Table of contents

| 1 | Components | | |
|---|------------|--|----|
| | 1.1 | Contents and storage | 5 |
| | 1.2 | Additional material to be supplied by user | 5 |
| 2 | Proc | duct description | 7 |
| | 2.1 | The basic principle | 7 |
| | 2.2 | Specifications | 8 |
| | 2.3 | General information | 10 |
| | 2.4 | Compatibility of reagents | 11 |
| 3 | Safe | aty instructions | 13 |
| | 3.1 | Risk and safety phrases | 13 |
| | 3.2 | GHS classification | 14 |
| 4 | Puri | fication of polyhistidine-tagged proteins under native conditions | 15 |
| | 4.1 | Preparation of buffers for purification under native conditions | 15 |
| | 4.2 | Preparation of cleared E. coli lysates under native conditions | 16 |
| | 4.3 | Batch purification of polyhistidine-tagged proteins under native conditions | 18 |
| | 4.4 | Semi-batch purification of polyhistidine-tagged proteins under native conditions | 20 |
| | 4.5 | Spin column purification of polyhistidine-tagged proteins under native conditions | 22 |
| | 4.6 | Gravity-flow purification of polyhistidine-tagged proteins under native conditions | 24 |
| | 4.7 | FPLC [™] purification of polyhistidine-tagged proteins under native conditions (self-packed columns) | 26 |
| | 4.8 | FPLC [™] purification of polyhistidine-tagged proteins under native conditions using Protino [®] Ni-NTA Columns 1 mL/5 mL | 28 |
| 5 | Puri | fication of polyhistidine-tagged proteins under denaturing conditions | 30 |
| | 5.1 | Preparation of buffers for purification under denaturing conditions | 30 |
| | 5.2 | Protein extract preparation under denaturing conditions | 32 |
| | 5.3 | Purification under denaturing conditions using Protino [®] Ni-NTA | 33 |
| 6 | Clea | aning, recharging, and storage | 34 |

| 7 | Арр | Appendix | |
|---|-----|------------------------------------|----|
| | 7.1 | Troubleshooting | 35 |
| | 7.2 | Ordering information | 37 |
| | 7.3 | Product use restriction / warranty | 38 |
| | | | |

1 Components

1.1 Contents and storage

| | Protino [®] Ni-NTA Agarose | | | |
|-------------------------------------|-------------------------------------|------------|------------|--|
| REF | 745400.25 | 745400.100 | 745400.500 | |
| Protino [®] Ni-NTA Agarose | 25 mL | 100 mL | 500 mL | |
| User manual | 1 | 1 | 1 | |

| | Protino [®] Ni-NTA Columns | | | |
|-------------------------------------|-------------------------------------|----------|----------|--|
| | 1 mL | 5 ו | mL | |
| REF | 745410.5 | 745415.1 | 745415.5 | |
| Protino [®] Ni-NTA Columns | 5 x 1 mL | 1 x 5 mL | 5 x 5 mL | |
| User manual | 1 | 1 | 1 | |

Shipping and storage

The product is shipped at ambient temperature. Upon receipt Protino[®] Ni-NTA Agarose products should be **stored at 4–8** °C and are stable up to 1 year. Do not freeze.

1.2 Additional material to be supplied by user

- For the purification under native conditions prepare the following buffers: NPI-10, NPI-20, NPI-250 (see section 4.1)
- For the purification under denaturing conditions prepare the following buffers: NPI-10, DNPI-10, DNPI-20, DNPI-250 (see section 5.1)
- Lysozyme
- Appropriate columns, centrifuge tubes, etc.
- Appropriate centrifuge

- Liquid chromatography system (MPLC, FPLC[™], ÄKTAdesign[™], etc.), peristaltic pump, or syringe
- If necessary, appropriate adaptors for connecting the Protino[®] Ni-NTA Columns to the system of choice. Protino[®] Ni-NTA Columns are equipped with 10–32 (1/16") inlet and outlet ports. With these ports the columns can easily be connected to standard MPLC/FPLC[™] systems (e.g., ÄKTAdesign[™]). Five adaptor sets are available for connecting the columns to other systems or for using them with a syringe.

Table 1: Adaptor sets

| System | Adaptor needed | Adaptor Set |
|---|--|---|
| Standard FPLC™ system (e.g., ÄKTAdesign™) | None | None |
| FPLC [™] system, first generation (Pharmacia) | 1 x M6 female to 10–32 male 1 x 10–32 female to M6 male | Protino [®] M6 Adaptor Set, REF 745260 |
| MPLC system (e.g., BioLogic [™] , BIO-RAD) | 1 x 1/4" 28 female to 10–32 male 1 x 10–32 female to 1/4" 28 female | Protino [®] 1/4-28 Adaptor Set, REF 745261 |
| MPLC system (e.g., Profinia™, BIO-RAD) | 1 x Luer female to 10–32 male 1 x 10–32 female to Luer male | Protino [®] Luer Adaptor Set, REF 745264 |
| Peristaltic pump | 1 x 1/16" ID tubing to 10–32 male | Protino [®] Inlet PP Adaptor Set, REF 745263 |
| Syringe | 1 x Luer female to 10–32 male | Protino [®] Inlet Luer Adaptor, REF 745262 |

2 Product description

2.1 The basic principle

Protino[®] Ni-NTA products enable fast and convenient purification of recombinant polyhistidine-tagged proteins by immobilized metal ion affinity chromatography (IMAC). Proteins from any expression system can be purified under native or denaturing conditions. Binding of protein is based on the interaction between the polyhistidine tag of the recombinant protein and immobilized Ni²⁺ ions.

Protino[®] Ni-NTA Agarose consists of the chelating ligand nitrilotriacetic acid (NTA) immobilized on 6 % cross-linked agarose beads that are suitable for batch binding, gravity flow, and FPLC[™] columns. The resin is precharged with Ni²⁺ ions and therefore ready to use.

Protino[®] Ni-NTA Agarose uses NTA which represents the most commonly used chelating ligand in IMAC. NTA is a tetradentate chelator which occupies four out of the six binding sites in the coordination sphere of the Ni²⁺ ion. The remaining two coordination sites are usually occupied by water molecules and can be exchanged with histidine residues of the recombinant protein (Figure 1). This formation of coordination sites has turned out to be optimal for purification of polyhistidine-tagged proteins: two available binding sites in the coordination sphere of the Ni²⁺ ion enable tight but reversible selective protein interactions. Chelation of Ni²⁺ ions by NTA through four coordination positions minimizes metal leaching during purification and increases specificity for polyhistidine-tagged proteins.

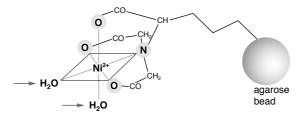


Figure 1: Protino[®] Ni-NTA Agarose – structure of NTA in complex with Ni²⁺

2.2 Specifications

| Table 2: Specifications Protino [®] Ni-NTA Agarose | | | | |
|---|--|--|--|--|
| Application | Batch binding Gravity-flow column chromatography MPLC/FPLC[™] | | | |
| Form | 50% aqueous suspension containing 30 vol% ethanol, precharged with Ni ²⁺ | | | |
| Support | Cross-linked 6 % beaded agarose | | | |
| Ligand | Nitrilotriacetic acid (NTA) | | | |
| Bead size | 45–165 μm | | | |
| Binding capacity ¹ | Up to 50 mg/mL settled agarose | | | |
| Recommended flow rate | <u>1 mL bed volume (column with 6.6 mm inner diameter)</u> 1.0 mL/min | | | |
| | <u>10 mL bed volume (column with 16 mm inner</u> diameter) 5.0 mL/min | | | |
| Maximum linear flow rate ² | 300 cm/h | | | |
| Storage temperature | 4–8 °C (do not freeze) | | | |

¹ Binding capacity will vary for each polyhistidine-tagged protein.

² High flow rates may reduce binding capacity.

| Table 3: Specifications Protino | [®] Ni-NTA Columns 1 mL/5 mL |
|---------------------------------|---------------------------------------|
|---------------------------------|---------------------------------------|

| Column bed volume | 1 mL | 5 mL | |
|--|---|---|--|
| System compatibility | Automated liquid chr (MPLC, FPLC[™], ÄK Peristaltic pump Syringe | omatography systems TAdesign™, etc.) | |
| Column dimensions | 0.7 cm inner diameter1.6 cm inner diameterx 2.5 cm heightx 2.5 cm height | | |
| Column body material | Polypropylene | | |
| Column ports | Inlet 10–32 (1/16") female Outlet 10–32 (1/16") male | | |
| Support | Support Cross-linked 6 % beaded agarose | | |
| Ligand | Nitrilotriacetic acid (NTA) | | |
| Bead size | 45–165 μm | | |
| Binding capacity ¹ | Up to 50 mg Up to 250 mg | | |
| Maximum back pressure | 3 bar (0.3 MPa) | | |
| Recommended flow rates | 1 mL/min | 5 mL/min | |
| Maximum recommended flow rate ² | ended 4 mL/min 10 mL/min | | |
| Storage temperature | 4–8 °C (do not freeze) | | |
| Storage solution | 30 vol % ethanol | | |

¹ Binding capacity will vary for each polyhistidine-tagged protein.

² High flow rates may reduce binding capacity.

2.3 General information

Binding capacity

- The binding capacity of Protino[®] Ni-NTA Agarose strongly depends on the characteristics of the polyhistidine-tagged protein (e.g., amino acid composition, molecular weight, 3-D structure, oligomerization properties). Furthermore, the absolute yield also depends on the total amount and concentration of the target protein in the sample which in turn directly correlate with the expression level and the cell density of the expression culture. Therefore binding capacity will vary for each polyhistidine-tagged protein.
- For best results determine the binding behaviour of any polyhistidine-tagged protein prior to attempting large-scale purification.
- A maximum capacity of up to 50 mg/mL was determined for the monomeric green fluorescent protein (6xHis-GFPuv, ~ 32 kDa) expressed in *E. coli*.

Solubility of the recombinant protein

- Protein yield is also dependent on solubility of the recombinant protein. If proteins are expressed in *E. coli*, ideally the target proteins remain soluble in the cytoplasm. However, especially proteins that are highly expressed accumulate in insoluble aggregates which are called inclusion bodies.
- For solubilization of inclusion bodies, buffers containing large amounts of denaturants are used.
- This manual includes instructions for isolation of soluble proteins (purification under native conditions, see section 4) as well as insoluble proteins from inclusion bodies (purification under denaturing conditions, see section 5).

Improving purity

- Sometimes optimization of purification procedures is necessary to increase purity.
- Usually lysis/equilibration buffers and the wash buffer contain 10 mM and 20 mM imidazole, respectively, to suppress binding of contaminating proteins. To improve specificity increase imidazole concentration.
- In addition, for more stringent binding and washing conditions the pH may be reduced from pH 8 closer to pH 7 (e.g., pH 7.4) in all buffers

Additives

• Avoid high concentration of additives that interact with nickel ions and thus reduce capacity (e.g., chelating agents (EDTA) or reducing agents (DTT, mercaptoethanol)), see compatibility of reagents (section 2.4).

2.4 Compatibility of reagents

| Table 4: Reagent compatibility chart | | | | | |
|--------------------------------------|---|---|--|--|--|
| Reagent | Effect | Comments | | | |
| Sodium phosphate | Used in buffers in order to buffer the solutions at pH 8 | 50 mM is recommended; the pH of any buffer should be adjusted to 8, although in some cases a pH between 7 and 8 can be used | | | |
| Tris, HEPES, MOPS | Coordinates with Ni ²⁺ ions, causing a decrease in capacity | Up to 100 mM may be used, sodium phosphate buffer is recommended | | | |
| Sodium Chloride | Prevents ionic interactions and therefore unspecific binding | Up to 2 M can be used, at least 0.3 M should be used | | | |
| Imidazole | Binds to immobilized Ni ²⁺ ions and competes with the polyhistidine-tagged proteins | Is used at low concentration to reduce non specific binding (20 mM) and to elute the target protein (> 100 mM) | | | |
| Urea | Solubilizes protein | Use 8 M for purification under denaturing conditions | | | |
| GuHCI | Solubilizes protein | Up to 6 M can be used | | | |
| ß-mercaptoethanol | Prevents formation of disulfide bonds; can reduce Ni ²⁺ ions at higher concentrations | Up to 20 mM in samples has been used successfully in some cases | | | |
| DTT, DTE | Can reduce Ni ²⁺ ions at higher concentrations | Up to 10 mM in samples has been used successfully in some cases | | | |
| Glutathione reduced | Can reduce Ni ²⁺ ions at higher concentrations | Up to 30 mM in samples has been used successfully in some cases | | | |
| Glycerol | Prevents hydrophobic interactions between proteins | Up to 50 % can be used | | | |
| EDTA | Coordinates with Ni ²⁺ ions, causing a decrease in capacity at higher concentrations | Not recommended, but up to 1 mM in samples has been used successfully in some cases | | | |

| Table 4: Reagent compatibility chart | | | | | |
|---|--|---|--|--|--|
| Ethanol | Prevents hydrophobic interactions between proteins | Up to 20 % can be used; ethanol may precipitate proteins, causing low flow rates and column clogging | | | |
| SDS | Interacts with Ni ²⁺ ions, causing a decrease in capacity | Not recommended, but up to 0.3 % in samples has been used successfully in some cases | | | |
| Nonionic detergents: Triton, Tween, etc. | Removes background proteins | Up to 2 % can be used | | | |

3 Safety instructions

The following components of the **Protino[®] Ni-NTA** products contain hazardous contents.

Wear gloves and goggles and follow the safety instructions given in this section.

3.1 Risk and safety phrases

| Component | Hazard contents | Hazard symbol | Risk phrases | Safety phrases |
|--|--|------------------------|-----------------|-------------------|
| Inhalt | Gefahrstoff | Gefahrstoff- symbol | R-Sätze | S-Sätze |
| Protino [®] Ni-NTA Agarose; Protino [®] Ni-NTA Columns 1 mL/5 mL | Ethanol 5–20 % + Ni in nickel(II) complexes 0.01–0.1 % Ethanol 5–20 % + Ni in Nickel(II)-Komplexen 0,01–0,1 % | Xi Xi | 10-43 | 37 |

Risk phrases

- R 10 Flammable. Entzündlich.
- R 43 May cause sensitisation by skin contact. Sensibilisierung durch Hautkontakt möglich.

Safety phrases

S 37 Wear suitable gloves. Geeignete Schutzhandschuhe tragen.

3.2 GHS classification

Only harmful features need not be labeled with H and P phrases until 125 mL or 125 g. *Mindergefährliche Eigenschaften müssen bis 125 mL oder 125 g nicht mit H- und P-Sätzen gekennzeichnet werden.*

| Component | Hazard contents | GHS symbol | Hazard phrases | Precaution phrases |
|--|-----------------|--------------------|-------------------|--|
| Inhalt | Gefahrstoff | GHS Symbol | H-Sätze | P-Sätze |
| Protino [®] Ni-NTA Agarose; Protino [®] Ni-NTA Columns 1 mL/5 mL | | Warning Achtung | 226, 317 | 261, 272, 280, 302+352, 333+313, 363 |

Hazard phrases

- H 226 Flammable liquid and vapour. Flüssigkeit und Dampf entzündbar.
- H 317 May cause an allergic skin reaction. Kann allergische Hautreaktionen verursachen.

Precaution phrases

| P 261 | Avoid breathing dust. Einatmen von Staub vermeiden. |
|-----------|---|
| P 272 | Contaminated work clothing should not be allowed out of the workplace. Kontaminierte Arbeitskleidung nicht außerhalb des Arbeitsplatzes tragen. |
| P 280 | Wear protective gloves / eye protection. Schutzhandschuhe / Augenschutz tragen. |
| P 302+352 | IF ON SKIN: Wash with plenty of soap and water. BEI KONTAKT MIT DER HAUT: Mit viel Wasser und Seife waschen. |
| P 333+313 | If skin irritation occurs: Get medical advice / attention. Bei Hautreizung: Ärztlichen Rat einholen / ärztliche Hilfe hinzuziehen. |
| P 363 | Wash contaminated clothing before reuse Kontaminierte Kleidung vor erneutem Tragen waschen. |

For further information please see Material Safety Data Sheets (www.mn-net.com). Weiterführende Informationen finden Sie in den Sicherheitsdatenblättern (www.mn-net.com).

A The symbol shown on labels refers to the precaution phrases of this section.

Das auf Etiketten dargestellte Symbol weist auf die P-Sätzen dieses Kapitels hin.

4 Purification of polyhistidine-tagged proteins under native conditions

4.1 Preparation of buffers for purification under native conditions

| NPI-10 / lysis & equilibra | tion buffer (1 liter): | |
|--|--|-------------------------------|
| $50 \text{ mM NaH}_2\text{PO}_4$ | 7.80 g NaH ₂ PO ₄ • 2 H ₂ O | M _r = 156.01 g/mol |
| 300 mM NaCl | 17.54 g NaCl | M _r = 58.44 g/mol |
| 10 mM imidazole | 0.68 g imidazole | M _r = 68.08 g/mol |
| Adjust pH to 8.0 using | NaOH | |
| NPI-20 / wash buffer (1 lit | er): | |
| 50 mM NaH ₂ PO ₄ | 7.80 g NaH ₂ PO ₄ • 2 H ₂ O | M _r = 156.01 g/mol |
| 300 mM NaCl | 17.54 g NaCl | M _r = 58.44 g/mol |
| 20 mM imidazole | 1.36 g imidazole | M _r = 68.08 g/mol |
| Adjust pH to 8.0 using | NaOH | |
| NPI-250 / elution buffer (* | l liter): | |
| 50 mM NaH ₂ PO ₄ | 7.80 g NaH ₂ PO ₄ • 2 H ₂ O | M _r = 156.01 g/mol |
| 300 mM NaCl | 17.54 g NaCl | M _r = 58.44 g/mol |
| 250 mM imidazole | 17.00 g imidazole | M _r = 68.08 g/mol |
| Adjust pH to 8.0 using | NaOH | |

4.2 Preparation of cleared *E. coli* lysates under native conditions

Cultivate and harvest cells

- Harvest cells from an *E. coli* expression culture by centrifugation at 4,500–6,000 x g for 15 min at 4 °C. Remove supernatant.
- To wash the cells resuspend in NPI-10 and centrifuge again. Remove supernatant.
- Cell pellets may be stored at -20 °C or -80 °C until needed.

Resuspend bacteria cells

- Thaw the cell pellet from an *E. coli* expression culture on ice (if frozen).
- Resuspend 1 g of pelleted, wet cells in 2–5 mL NPI-10. Pipette up and down, or stir until complete resuspension without visible cell aggregates. Perform this step on ice.

Lyse cells

- Add lysozyme to a final concentration of 1 mg/mL.
- Stir the solution on ice for 30 min.
- Sonicate the suspension on ice according to the instructions provided by the manufacturer (e.g., use 10 x 15 s bursts with a 15 s cooling period between each burst).
- Carefully check samples' appearance after sonication. If the lysate is still viscous from incomplete fragmentation of DNA, add 5 µg/mL DNase I and stir on ice for 15 min.

Clarify lysate

- Centrifuge the crude lysate at 10,000 x g for 30 min at 4 °C to remove cellular debris.
- Carefully transfer the supernatant to a clean tube without disturbing the pellet. If the supernatant is not clear, centrifuge a second time or filter through a 0.45 µm membrane (e.g., cellulose acetate).
- Store supernatant on ice.

| Proceed to section | 4.3 for batch purification |
|--------------------|---|
| | 4.4 for semi-batch purification |
| | 4.5 for Spin column purification |
| | 4.6 for gravity-flow purification |
| | 4.7 for FPLC [™] purification using self-packed columns or |
| | 4.8 for FPLC [™] purification using Protino [®] Ni-NTA Columns. |

4.3 Batch purification of polyhistidine-tagged proteins under native conditions

1 Equilibration

- Resuspend Protino[®] Ni-NTA Agarose by mixing thoroughly to achieve a homogeneous suspension.
- Immediately transfer an appropriate amount of suspension to an appropriate tube.

Pipette 2 mL of the original 50 % suspension per 1 mL of bed volume required.

- Sediment the gel by centrifugation at 500 x *g* for 5 minutes. Carefully decant the supernatant (storage solution) and discard it.
- Add 10 bed volumes of NPI-10 to equilibrate the gel. Invert to mix.
- Sediment the gel by centrifugation at 500 x *g* for 5 minutes. Carefully decant the supernatant and discard it.

2 Batch binding

- Add the clarified *E. coli* lysate or protein extract to the equilibrated gel.
- Mix the suspension gently for 30-60 min.
- Sediment the gel by centrifugation at 500 x *g* for 5 minutes. Carefully decant the supernatant and discard it.

3 Washing

- Wash the gel by adding 10 bed volumes of NPI-20. Invert to mix.
- Sediment the gel by centrifugation at 500 x *g* for 5 minutes. Carefully decant the supernatant and discard it.
- Repeat the washing step (total wash 2 x 10 bed volumes of NPI-20).

4 Elution

- Add 1 bed volume of NPI-250 to the sedimented gel.
- Mix the suspension gently for 2 min at room temperature to liberate the polyhistidine-tagged protein from the gel.
- Sediment the gel by centrifugation at 500 x *g* for 5 minutes. Carefully decant or pipette the supernatant in a new tube and store eluted protein on ice.

- Repeat the elution step 5 times.
- Analyze fractions for the presence of the target protein. To determine the protein concentration use a Bradford protein assay (quick and easy) or measure the absorbance at 280 nm. Note that imidazole will also absorb at 280 nm.
- Pool fractions containing the majority of the eluted polyhistidine-tagged protein and remove an aliquot for SDS-PAGE analysis.
- Store protein at -70 °C or -20 °C. Note that many proteins irreversibly precipitate out of solution in the presence of imidazole. In this case remove imidazole prior to freezing.

For proper storage it is recommended to remove the imidazole by ultrafiltration or dialysis.

4.4 Semi-batch purification of polyhistidine-tagged proteins under native conditions

1 Equilibration

- Resuspend Protino[®] Ni-NTA Agarose by mixing thoroughly to achieve a homogeneous suspension
- Immediately transfer an appropriate amount of suspension to an appropriate chromatography column.

Pipette 2 mL of the original 50% suspension per 1 mL of bed volume required.

- Allow the column to drain by gravity.
- Add 10 bed volumes of NPI-10 to equilibrate the gel.
- Allow the column to drain by gravity.

2 Batch binding

- Close column outlet with cap.
- Add the clarified *E. coli* lysate or protein extract to the equilibrated gel.
- · Close column inlet with a cap.
- Mix the suspension gently for 30–60 min by slowly inverting the column.
- Install the column in a vertical position.
- · Remove bottom and top caps.
- Allow the column to drain by gravity

3 Washing

- · Wash the column with 10 bed volumes of NPI-20
- · Allow the column to drain by gravity.
- Repeat the washing step (total wash 2 x 10 bed volumes of NPI-20).

4 Elution

- Add 5–10 bed volumes of NPI-250 to the gel.
- Allow the column to drain by gravity and collect the eluate in fractions.
- · Store eluted protein on ice.

- Analyze fractions for the presence of the target protein. To determine the protein concentration use a Bradford protein assay (quick and easy) or measure the absorbance at 280 nm. Note that imidazole will also absorb at 280 nm.
- Pool fractions containing the majority of the eluted polyhistidine-tagged protein and remove an aliquot for SDS-PAGE analysis.
- Store protein at -70 °C or -20 °C. Note that many proteins irreversibly precipitate out of solution in the presence of imidazole. In this case remove imidazole prior to freezing.

For proper storage it is recommended to remove the imidazole by ultrafiltration or dialysis.

4.5 Spin column purification of polyhistidine-tagged proteins under native conditions

1 Equilibration

- Resuspend Protino[®] Ni-NTA Agarose by mixing thoroughly to achieve a homogeneous 50 % suspension.
- Immediately transfer 50 μL of the suspension to a Receiver Column placed in a collecting tube.
- Centrifuge at 500 x g for 30 s.
- To equilibrate the gel add 500 μL of NPI-10.
- Centrifuge at 500 x g for 30 s. Discard flow-through.

2 Binding

- Close column outlet with cap.
- Add up to 700 µL of clarified *E. coli* lysate to the gel and close the lid.
- Mix the suspension gently for 1 h at room temperature, for example, using an eppendorf thermomixer.
- Remove bottom cap and place Receiver Column in a collecting tube.
- Centrifuge at 500 x g for 30 s. Discard flow-through.

3 Washing

- To wash the gel add 500 μL of NPI-20.
- Centrifuge at 500 x g for 30 s. Discard flow-through.
- Repeat the washing step twice (total wash 3 x 500 μL of NPI-20). Discard flow-through between washing steps.

4 Elution

• Close column outlet with cap. Add 125–250 μL of NPI-250 and close the lid.

<u>Note</u>: The amount of elution buffer required depends on the amount of loaded protein and may vary among fusion proteins. Volumes may have to be adjusted accordingly.

- Mix the suspension gently for 2 min at room temperature to liberate polyhistidine-tagged proteins from the gel.
- Remove bottom cap and place Receiver Column in a 1.5 or 2 mL microcentrifuge tube.
- Centrifuge at 500 x g for 30 s.
- Repeat the elution step 2 times.
- Store eluted protein on ice.
- Analyze fractions for the presence of the target protein. To determine the protein concentration use a Bradford protein assay (quick and easy) or measure the absorbance at 280 nm. Note that imidazole will also absorb at 280 nm.
- Pool fractions containing the majority of the eluted polyhistidine-tagged protein and remove an aliquot for SDS-PAGE analysis.
- Store protein at -70 °C or -20 °C. Note that many proteins irreversibly precipitate out of solution in the presence of imidazole. In this case remove imidazole prior to freezing.
- For proper storage conditions or downstream application imidazole can be removed by gelfiltration, ultrafiltration or dialysis.

4.6 Gravity-flow purification of polyhistidine-tagged proteins under native conditions

1 Equilibration

- Resuspend Protino[®] Ni-NTA Agarose by mixing thoroughly to achieve a homogeneous suspension.
- Immediately transfer an appropriate amount of suspension to an appropriate chromatography column, which allows slow flow rates of 0.5–1 mL/min.

Pipette 2 mL of the original 50 % suspension per 1 mL of bed volume required.

- Allow the column to drain by gravity.
- Add 10 bed volumes of NPI-10 to equilibrate the gel.
- Allow the column to drain by gravity.

2 Binding

- Add the clarified *E. coli* lysate or protein extract to the equilibrated gel.
- Allow the column to drain by gravity using a flow rate of 0.5–1 mL/min.

<u>Note</u>: If the flow rate is too high polyhistidine-tagged proteins may not bind to the column efficiently. Reduce the flow rate or re-apply the flow-through to improve binding.

3 Washing

- Wash the gel by adding 10 bed volumes of NPI-20.
- Allow the column to drain by gravity.
- Repeat the washing step (total wash 2 x 10 bed volumes of NPI-20).

4 Elution

- Add 5–10 bed volumes of NPI-250 to the gel.
- Allow the column to drain by gravity and collect the eluate in fractions.
- · Store eluted protein on ice.
- Analyze fractions for the presence of the target protein. To determine the protein concentration use a Bradford protein assay (quick and easy) or measure the absorbance at 280 nm. Note that imidazole will also absorb at 280 nm.

- Pool fractions containing the majority of the eluted polyhistidine-tagged protein and remove an aliquot for SDS-PAGE analysis.
- Store protein at -70 °C or -20 °C. Note that many proteins irreversibly precipitate out of solution in the presence of imidazole. In this case remove imidazole prior to freezing.

For proper storage it is recommended to remove the imidazole by ultrafiltration or dialysis.

4.7 FPLC[™] purification of polyhistidine-tagged proteins under native conditions (self-packed columns)

Prepare buffers according to section 4.1. Filter buffers through a 0.45 μm filter before use. Clear samples by centrifugation and/or pass them through a 0.45 μm filter.

1 Preparing the chromatography system

- Purge the pump with deionized water. Assure that all air is displaced.
- Determine the bed volume of Protino[®] Ni-NTA Agarose required for your application. Choose an appropriate chromatography column (e.g., from Omnifit or GE Healthcare). If more than 50 % of the column volume is to be packed, equip the column with an extension to hold the complete volume of the agarose suspension.
- Eliminate air from outlet tubing and end piece of the column by injecting deionized water into outlet tubing. Close outlet of column. Leave ~ 1 cm of buffer above the support net or frit.
- Inject deionized water into the inlet tubing of the upper plunger to eliminate air. Place plunger into a beaker containing deionized water until use.

2 Column packing

- Resuspend Protino[®] Ni-NTA Agarose by mixing thoroughly to achieve a homogeneous 50% suspension. Immediately transfer the determined volume of suspension to an appropriate vacuum flask and de-gas.
- Pour the entire slurry into the column in one continuous motion along a glass rod held against the inner wall of the column.
- Carefully fill the remaining space with deionized water. Insert the upper plunger into the column without introducing air bubbles. Connect the inlet of the column to a pump.
- Open the column outlet and start the pump. Pass deionized water through the column at a packing flow rate of approximately 300 cm/h until height of gel bed becomes constant. Stop the pump and close the column outlet.
- Position the upper plunger on top of the column bed. Avoid to introduce air bubbles. Open the column outlet and start the pump at a flow rate of approximately 300 cm/h until the bed is stable. Re-position the plunger on the medium surface as necessary.

3 Column equilibration

- Purge the pump with NPI-10.
- Equilibrate the column with 5–10 bed volumes of NPI-10 until the baseline at 280 nm is stable.

4 Binding

- Load the clarified *E. coli* lysate or protein extract onto the column.
- Collect flow-through and analyze (e.g., by SDS-PAGE) to verify that the polyhistidine-tagged protein has bound. If the fusion protein is found in early fractions of the flow-through, the flow rate should be decreased. If the fusion protein is absent in early fractions and does appear in late fractions of the flow-through the column capacity has been exceeded. In this case protein load should be reduced or bed volume should be increased.

5 Washing

• Wash the column with 10–20 bed volumes of NPI-20 or until the baseline at 280 nm is stable.

6 Elution

- Elute the polyhistidine-tagged protein with 5–10 bed volumes of NPI-250 and collect fractions.
- Store eluted protein on ice.
- Analyze fractions for the presence of the target protein. If a 280 nm absorbance flow monitor is not available, determine the protein concentration using a Bradford protein assay (quick and easy) or measure the absorbance at 280 nm. Note that imidazole will also absorb at 280 nm.
- Pool fractions containing the majority of the eluted polyhistidine-tagged protein and remove an aliquot for SDS-PAGE analysis.
- Store protein at -70 °C or -20 °C. Note that many proteins irreversibly precipitate out of solution in the presence of imidazole. In this case remove imidazole prior to freezing.
- For proper storage it is recommended to remove the imidazole by ultrafiltration or dialysis.

4.8 FPLC[™] purification of polyhistidine-tagged proteins under native conditions using Protino[®] Ni-NTA Columns 1 mL/5 mL

Protino[®] Ni-NTA Columns can be operated with liquid chromatography systems (such as ÄKTAdesign[™] systems) via standard 10–32 fittings without additional connectors.

Prepare buffers according to section 4.1. Filter buffers through a 0.45 μm filter before use. Clear samples by centrifugation and/or pass them through a 0.45 μm filter

Protino[®] Ni-NTA Columns




1 Connect column to the chromatography system

- Purge the pump with NPI-10. Assure that all air is displaced.
- Remove the snap-off end at the column outlet and save it for further use.
- Remove the upper plug from the column.
- Start the pump at a flow rate of approximately 0.3 mL/min.
- Fill the inlet port of the column with several drops of NPI-10 to remove air to form a positive meniscus.
- Insert the fitting "drop-to-drop" into the column port to avoid introducing air bubbles.

<u>Note</u>: The snap-off end can be reused as a stop plug for sealing the column outlet for storage.

2 Column equilibration

 Equilibrate the column with 5–10 column volumes of NPI-10 until the baseline at 280 nm is stable.

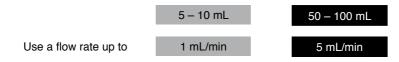

3 Binding

- Load the centrifuged or filtered sample onto the column.
- Use a flow rate up to
 1.0 mL/min
 5 mL/min

Protino[®] Ni-NTA Columns



 Collect flow-through and analyze (e.g., by SDS-PAGE) to verify that the polyhistidine-tagged protein has bound. If the fusion protein is found in early fractions of the flow-through, the flow rate should be decreased. If the fusion protein is absent in early fractions and does appear in late fractions of the flow-through, the column capacity has been exceeded. In this case protein load should be reduced or bed volume should be increased.

4 Washing

• Wash the column with 10–20 column volumes of NPI-20 or until the baseline at 280 nm is stable.



5 Elution

 Elute the polyhistidine-tagged protein with 5–10 column volumes of NPI-250 and collect fractions.



- Store eluted protein on ice.
- Analyze fractions for the presence of the target protein. If a 280 nm absorbance flow monitor is not available, y determine the protein concentration use a Bradford protein assay (quick and easy) or measure the absorbance at 280 nm. Note that imidazole will also absorb at 280 nm.
- Pool fractions containing the majority of the eluted polyhistidine-tagged protein and remove an aliquot for SDS-PAGE analysis.
- Store protein at -70 °C or -20 °C. Note that many proteins irreversibly precipitate out of solution in the presence of imidazole. In this case remove imidazole prior to freezing.
- For proper storage it is recommended to remove the imidazole by ultrafiltration or dialysis.

5 Purification of polyhistidine-tagged proteins under denaturing conditions

We recommend this protocol if expression leads to the formation of inclusion bodies. Cells are disrupted under native conditions using lysozyme together with sonication. After centrifugation the polyhistidine-tagged protein is extracted and solubilized from the pellet by using a denaturant (8 M urea). The extract obtained is clarified by centrifugation and applied to Protino[®] Ni-NTA Agarose or Protino[®] Ni-NTA Columns under denaturing conditions.

5.1 Preparation of buffers for purification under denaturing conditions

| NPI-10 (1 liter): | | |
|----------------------------------|--|-------------------------------|
| $50 \text{ mM NaH}_2\text{PO}_4$ | 7.80 g NaH ₂ PO ₄ • 2 H ₂ O | M _r = 156.01 g/mol |
| 300 mM NaCl | 17.54 g NaCl | M _r = 58.44 g/mol |
| 10 mM imidazole | 0.68 g imidazole | M _r = 68.08 g/mol |
| Adjust pH to 8.0 using N | NaOH | |
| DNPI-10 (1 liter): | | |
| $50 \text{ mM NaH}_2\text{PO}_4$ | 7.80 g NaH ₂ PO ₄ • 2 H ₂ O | M _r = 156.01 g/mol |
| 300 mM NaCl | 17.54 g NaCl | M _r = 58.44 g/mol |
| 10 mM imidazole | 0.68 g imidazole | M _r = 68.08 g/mol |
| 8 M urea | 480 g urea | M _r = 60.06 g/mol |
| Adjust pH to 8.0 using N | NaOH | |
| DNPI-20 (1 liter): | | |
| $50 \text{ mM NaH}_2\text{PO}_4$ | 7.80 g NaH ₂ PO ₄ • 2 H ₂ O | M _r = 156.01 g/mol |
| 300 mM NaCl | 17.54 g NaCl | M _r = 58.44 g/mol |
| 20 mM imidazole | 1.36 g imidazole | M _r = 68.08 g/mol |
| 8 M urea | 480 g urea | M _r = 60.06 g/mol |
| Adjust pH to 8.0 using N | NaOH | |

.....

| DNPI-250 (1 liter): | | | | | |
|---------------------|-------------|------------|--------------------------|---------|--------------|
| 50 mM | NaH_2PO_4 | 7.80 g N | $la_2HPO_4 \cdot 2 H_2O$ | $M_r =$ | 156.01 g/mol |
| 300 mM | NaCl | 17.54 g N | laCl | $M_r =$ | 58.44 g/mol |
| 250 mM | imidazole | 17.00 g ir | nidazole | $M_r =$ | 68.08 g/mol |
| 8 M | urea | 480 g u | rea | $M_r =$ | 60.06 g/mol |

Adjust pH to 8.0 using NaOH

5.2 Protein extract preparation under denaturing conditions

1 Isolation of inclusion bodies

- Thaw the cell pellet from an *E. coli* expression culture on ice (if frozen). Resuspend 1 g of pelleted, wet cells in 5 mL NPI-10 buffer (without denaturant) on ice. Pipette up and down, or use stirring until complete resuspension without visible cell aggregates.
- Add lysozyme to a final concentration of 1 mg/mL. Stir the solution on ice for 30 min.
- Sonicate the suspension on ice according to the instructions provided by the manufacturer (e.g., use 10 x 15 s bursts with a 15 s cooling period between each burst).
- Carefully check samples' appearance after sonication. If the lysate is still viscous from incomplete fragmentation of DNA, add 5 µg/mL DNase I and stir on ice for 15 min
- Centrifuge the crude lysate at 10,000 x g for 30 min at 4 °C to collect the inclusion bodies. Discard supernatant. Keep pellet on ice.

2 Solubilization of inclusion bodies

- To wash the inclusion bodies resuspend the pellet in 10 mL NPI-10 (without denaturant) per g wet cells.
- Centrifuge suspension at 10,000 x g for 30 min at 4 °C. Discard supernatant.
- Resuspend the pellet in 2.0 mL DNPI-10 per g wet cells to solubilize the inclusion bodies. Homogenization or sonication may be necessary to resuspend the pellet. Dissolve the inclusion bodies by stirring on ice for 60 min.
- Centrifuge at 10,000 x g for 30 min at 20 °C to remove any remaining insoluble material. Carefully transfer the supernatant to a clean tube without disturbing the pellet.
- f the supernatant is not clear centrifuge a second time or filter through a 0.45 µm membrane (e.g., celluloseacetate) to avoid clogging of the IMAC column with insoluble material.
- Save supernatant (solubilized protein).

5.3 Purification under denaturing conditions using Protino[®] Ni-NTA

Purification of polyhistidine-tagged proteins under denaturing conditions is similar to purification under native conditions except that the sample and buffers loaded on the column contain 8 M urea.

1 Proceed to section 4.3 for batch purification

4.4 for semi-batch purification

4.5 for Spin column purification

4.6 for gravity-flow purification

- 4.7 for FPLC[™] purification using self-packed columns or
- 4.8 for FPLC[™] purification using Protino[®] Ni-NTA Columns

with the following modifications - use:

- · Supernatant from 5.2 (solubilized protein) as sample or protein extract,
- DNPI-10 instead of NPI-10 (equilibration buffer),
- DNPI-20 instead of NPI-20 (wash buffer),
- DNPI-250 instead of NPI-250 (elution buffer).

6 Cleaning, recharging, and storage

Cleaning

After use Protino[®] Ni-NTA Agarose should be washed for 30 minutes with 0.5 M NaOH followed by equilibration. We recommend this cleaning procedure if you wish to purify the same His-tag protein. Sodium hydroxide effectively desorbs contaminants originating from the loaded sample, such as unspecifically bound proteins, precipitated proteins and lipoproteins.

 Wash Protino[®] Ni-NTA Agarose with 15 bed volumes of 0.5 M NaOH for 30 min to solubilize and desorb contaminants.

When using columns adjust the flow rate accordingly. For example, wash a Protino[®] Ni-NTA Columns 1 mL by using a flow rate of 0.5 mL/min for 30 min, corresponding to a total volume of 15 mL.

- Remove the NaOH solution by washing with 10 bed volumes of de-ionized water.
- If you are reusing the resin directly, wash with 10 bed volumes of NPI-10 to equilibrate the resin.
- For storage wash with 2 bed volumes of 30% ethanol. Resuspend the resin in 30% ethanol and store at 2–8 °C.

Recharging

Depending on the nature of the sample the cleaning procedure mentioned above may not be satisfactory. In cases, for example when the color of the resin changes (from light blue to white/grey (due to loss of nickel ions) or to brown (due to the reduction of nickel ions)) Protino[®] Ni-NTA Agarose can easily be stripped and recharged with nickel.

- Wash Protino® Ni-NTA Agarose with 10 bed volumes of de-ionized water.
- Strip of nickel ions by washing with 10 bed volumes of 100 mM EDTA, pH 8.
- · Wash resin with 10 bed volumes of de-ionized water.
- Charge resin with 2 bed volumes of 100 mM metal ion aqueous solution (e.g. NiSO₄ or NiCl₂).

Other metal ions may be used to increase specificity (e.g., Co²⁺ or or Zn²⁺).

- Wash resin with 10 bed volumes of de-ionized water to remove unbound metal ions.
- If you are reusing the resin, directly wash with 10 bed volumes of NPI-10 to equilibrate the resin.
- For storage wash with 2 bed volumes of 30% ethanol. Resuspend the resin in 30% ethanol and store at 2–8 °C.

7 Appendix

7.1 Troubleshooting

| Problem | Possible cause and suggestions | | |
|--|---|--|--|
| Column is clogged | Sample /lysate contains insoluble material If the sample is not clear use centrifugation or filtration (0.45 μm membrane) to avoid clogging of the IMAC column. Sample /lysate contains genomic DNA Lysate may remain viscous from incomplete shearing of genomic DNA after sonication. Add 5 μg/mL DNase I and incubate on ice for 10 min. | | |
| Protein does not bind to the resin | Problems with vector construction Ensure that protein and tag are in frame. Sometimes the position of the tag influences expression rate and solubility. Evaluate N- and C-terminally tagged variants of the protein. His-Tag is not accessible. Use denaturing conditions to purify the protein. Use a C-terminal Histag instead of a N-terminal tag or vice versa. Incorrect binding conditions Check composition and pH of all buffers. Ensure that all additives are compatible (see compatibility of reagents, 2.4) | | |
| Protein elutes with wash buffer | Incorrect buffer compositionCheck composition and pH of all buffers. | | |
| Protein does not elute | Elution conditions are too mild. Increase concentration of imidazole from 250 mM to 500 mM. Protein has precipitated Elute under denaturing conditions. | | |

| Insufficient wash Use larger volumes for washing step. Use NPI-50 for third washing step (containing 50 mL in Binding and wash conditions are too mild Use 10–20 mM imidazole in the binding and washing Contaminating proteins and target protein are linked too disulfide bonds Add up to 20 mM 2-mercaptoethanol to reduce disulfide Contaminating proteins are proteolytic products of target Perform cell lysis at 4 °C. Include protease inhibitors. Poly- histidine- tagged protein Contaminating host proteins have a better chance to b resin when only small amounts of target protein are preselysate. Very low amounts of polyhistidine-tagged protei Protein Contaminating host proteins have a better chance to b resin when only small amounts of target protein are preselysate. Very low amounts of polyhistidine-tagged protei Preduce the amount of Protino® Ni-NTA resin or incluation of sample. Expression is too low Increase expression level. Sometimes the position influences expression rate and solubility. Use a of Histag instead of a N-terminal tag or vice versa. | |
|---|---|
| influences expression rate and solubility. Use a | g and washing buffers. are linked together via reduce disulfide bonds. ducts of target protein ein er chance to bind to the rotein are present in the -tagged protein are not ng proteins effectively. A resin or increase the |
| Increase amount of starting cell material.Do not exceed recommended lysis volumes. | bility. Use a C-terminal ce versa. al. |

7.2 Ordering information

| Product | REF | Pack of |
|--|---------------------------------------|------------------------------|
| Protino [®] Ni-NTA Agarose | 745400.25 745400.100 745400.500 | 25 mL 100 mL 500 mL |
| Protino [®] Ni-NTA Columns 1 mL | 745410.5 | 5 columns |
| Protino [®] Ni-NTA Columns 5 mL | 745415.1 745415.5 | 1 column 5 columns |
| Protino [®] 96 Ni-NTA | 745425.1 745425.4 | 1 x 96 preps 4 x 96 preps |
| Protino [®] Columns 14 mL (empty gravity-flow columns) | 745250.10 | 10 columns |
| Protino [®] Columns 35 mL (empty gravity-flow columns) | 745255.10 | 10 columns |
| Protino [®] M6 Adaptor Set | 745260 | 1 |
| Protino [®] 1/4-28 Adaptor Set | 745261 | 1 |
| Protino [®] Luer Adaptor Set | 745264 | 1 |
| Protino [®] Inlet PP Adaptor Set | 745263 | 1 |
| Protino [®] Inlet Luer Adaptor | 745262 | 1 |

Visit *www.mn-net.com* for more detailed product information.

7.3 Product use restriction / warranty

Protino[®] products components are intended, developed, designed, and sold FOR RESEARCH PURPOSES ONLY, except, however, any other function of the product being expressly described in original MACHEREY-NAGEL product leaflets.

MACHEREY-NAGEL products are intended for GENERAL LABORATORY USE ONLY! MACHEREY-NAGEL products are suited for QUALIFIED PERSONNEL ONLY! MACHEREY-NAGEL products shall in any event only be used wearing adequate PROTECTIVE CLOTHING. For detailed information please refer to the respective Material Safety Data Sheet of the product! MACHEREY-NAGEL products shall exclusively be used in an ADEQUATE TEST ENVIRONMENT. MACHEREY-NAGEL does not assume any responsibility for damages due to improper application of our products in other fields of application. Application on the human body is STRICTLY FORBIDDEN. The respective user is liable for any and all damages resulting from such application.

DNA/RNA/PROTEIN purification products of MACHEREY-NAGEL are suitable for *IN VITRO*-USES ONLY!

ONLY MACHEREY-NAGEL products specially labeled as IVD are also suitable for *IN VITRO*-diagnostic use. Please pay attention to the package of the product. *IN VITRO*-diagnostic products are expressly marked as IVD on the packaging.

IF THERE IS NO IVD SIGN, THE PRODUCT SHALL NOT BE SUITABLE FOR IN VITRO-DIAGNOSTIC USE!

ALL OTHER PRODUCTS NOT LABELED AS IVD ARE NOT SUITED FOR ANY CLINICAL USE (INCLUDING, BUT NOT LIMITED TO DIAGNOSTIC, THERAPEUTIC AND/OR PROGNOSTIC USE).

No claim or representations is intended for its use to identify any specific organism or for clinical use (included, but not limited to diagnostic, prognostic, therapeutic, or blood banking). It is rather in the responsibility of the user or - in any case of resale of the products - in the responsibility of the reseller to inspect and assure the use of the DNA/RNA/protein purification products of MACHEREY-NAGEL for a well-defined and specific application.

MACHEREY-NAGEL shall only be responsible for the product specifications and the performance range of MN products according to the specifications of in-house quality control, product documentation and marketing material.

This MACHEREY-NAGEL product is shipped with documentation stating specifications and other technical information. MACHEREY-NAGEL warrants to meet the stated specifications. MACHEREY-NAGEL's sole obligation and the customer's sole remedy is limited to replacement of products free of charge in the event products fail to perform as warranted. Supplementary reference is made to the general business terms and conditions of MACHEREY-NAGEL, which are printed on the price list. Please contact us if you wish to get an extra copy.

There is no warranty for and MACHEREY-NAGEL is not liable for damages or defects arising in shipping and handling (transport insurance for customers excluded), or out of accident or improper or abnormal use of this product; defects in products or components not manufactured by MACHEREY-NAGEL, or damages resulting from such non-MACHEREY-NAGEL components or products.

MACHEREY-NAGEL makes no other warranty of any kind whatsoever, and SPECIFICALLY DISCLAIMS AND EXCLUDES ALL OTHER WARRANTIES OF ANY KIND OR NATURE WHATSOEVER, DIRECTLY OR INDIRECTLY, EXPRESS OR IMPLIED, INCLUDING, WITHOUT LIMITATION, AS TO THE SUITABILITY, REPRODUCTIVITY, DURABILITY, FITNESS FOR A PARTICULAR PURPOSE OR USE, MERCHANTABILITY, CONDITION, OR ANY OTHER MATTER WITH RESPECT TO MACHEREY-NAGEL PRODUCTS.

In no event shall MACHEREY-NAGEL be liable for claims for any other damages, whether direct, indirect, incidental, compensatory, foreseeable, consequential, or special (including but not limited to loss of use, revenue or profit), whether based upon warranty, contract, tort (including negligence) or strict liability arising in connection with the sale or the failure of MACHEREY-NAGEL products to perform in accordance with the stated specifications. This warranty is exclusive and MACHEREY-NAGEL makes no other warranty expressed or implied.

The warranty provided herein and the data, specifications and descriptions of this MACHEREY-NAGEL product appearing in MACHEREY-NAGEL published catalogues and product literature are MACHEREY-NAGEL's sole representations concerning the product and warranty. No other statements or representations, written or oral, by MACHEREY-NAGEL's employees, agent or representatives, except written statements signed by a duly authorized officer of MACHEREY-NAGEL are authorized; they should not be relied upon by the customer and are not a part of the contract of sale or of this warranty.

Product claims are subject to change. Therefore please contact our Technical Service Team for the most up-to-date information on MACHEREY-NAGEL products. You may also contact your local distributor for general scientific information. Applications mentioned in MACHEREY-NAGEL literature are provided for informational purposes only. MACHEREY-NAGEL does not warrant that all applications have been tested in MACHEREY-NAGEL laboratories using MACHEREY-NAGEL products. MACHEREY-NAGEL does not warrant the correctness of any of those applications.

Last updated: 07/2010, Rev. 03

Please contact: MACHEREY-NAGEL GmbH & Co. KG Tel.: +49 24 21 969-270 tech-bio@mn-net.com

Trademarks:

ÄKTAdesig[™] and FPLC[™] are trademarks of GE Healthcare companies BioLogic[™] and Profinia[™] are trademarks of BIO-RAD Laboratories, Inc. Protino[®] is a registered trademark of MACHEREY-NAGEL

All used names and denotations can be brands, trademarks, or registered labels of their respective owner – also if they are not special denotation. To mention products and brands is only a kind of information (i.e., it does not offend against trademarks and brands and can not be seen as a kind of recommendation or assessment). Regarding these products or services we can not grant any guarantees regarding selection, efficiency, or operation.